
http://www.adsequip.com/

N4 Dental Stool
User Manual

P/N:8500900 Version：B Date：2024-10-31

1 Compositions

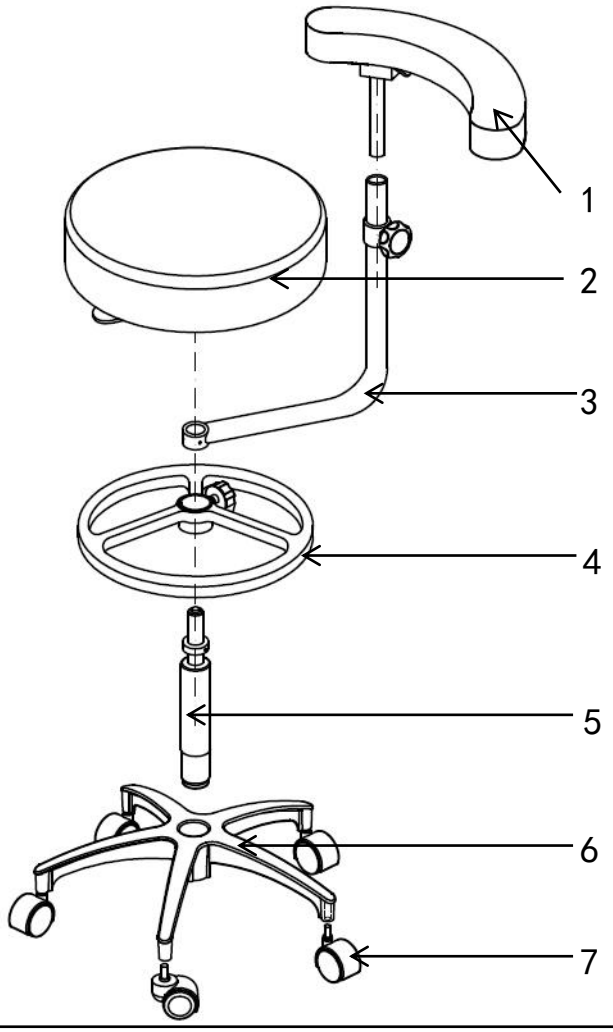


Item	Component
1	Armrest
2	Stool Seat
3	Mounting Arm Upright
4	Foot Ring
5	Height Adjustment Cylinder
6	Stool Base
7	Caster

Weight Load: 440 lbs


2 Installation

1. Install the casters:
Insert the casters into the stool base until fully seated.

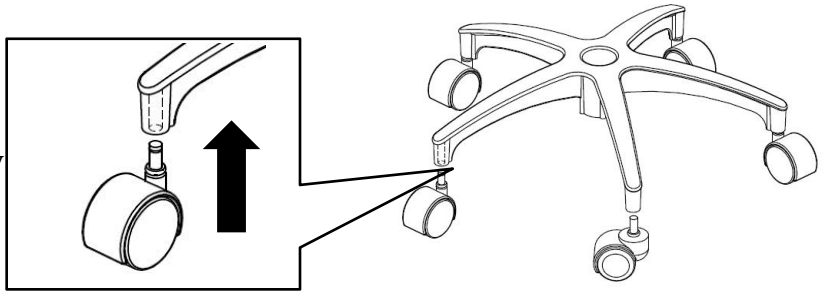

2. Install the height adjustable cylinder:


Insert the cylinder into the hole of the stool base until fully seated.

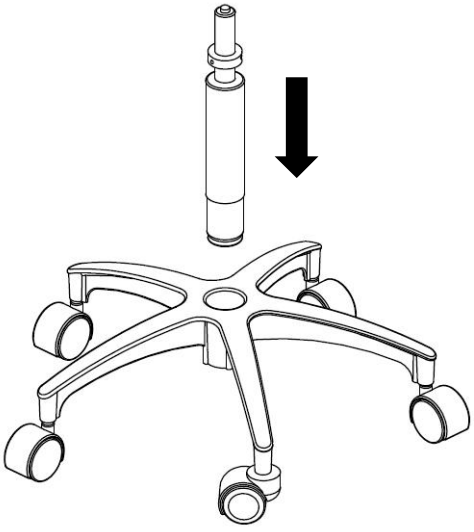

3 Install the foot ring:


Put the foot ring into the height adjustment cylinder.

（**Note:** If the mounting handwheel affects the installation of the foot ring, unscrew it.）

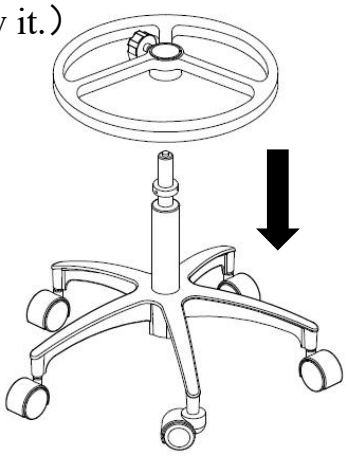


4. Install the armrest mounting arm upright:

Put the mounting arm upright into the height adjustment cylinder.

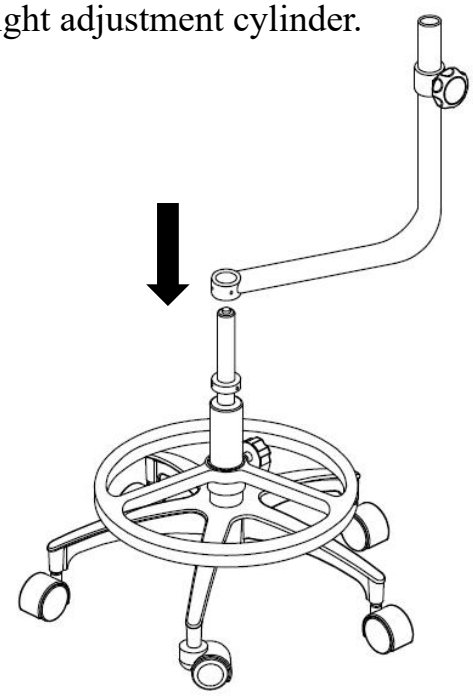


5. Install the stool seat:

Put the stool seat into the height adjustment cylinder.

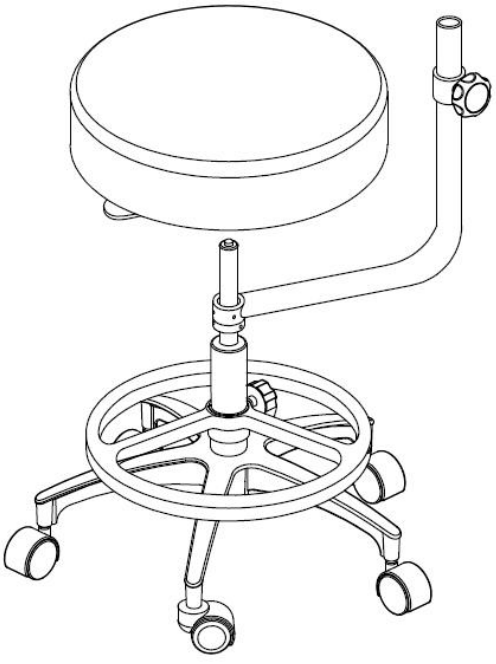


6. Install the armrest

Put the mounting arm upright of armrest into the hole and tighten the handwheel until fully seated.

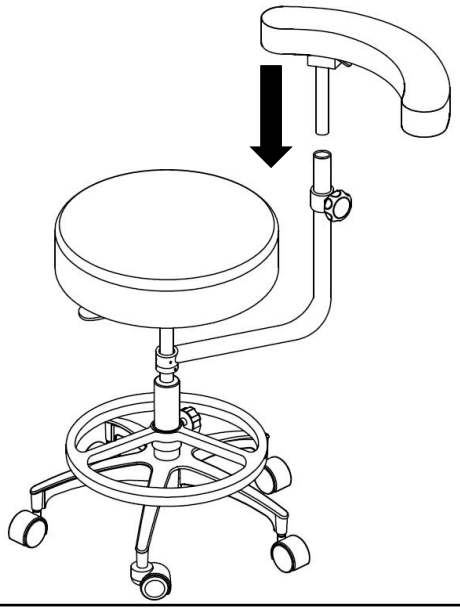


3 Adjustments

3.1 Adjust the Armrest Mounting Arm Upright

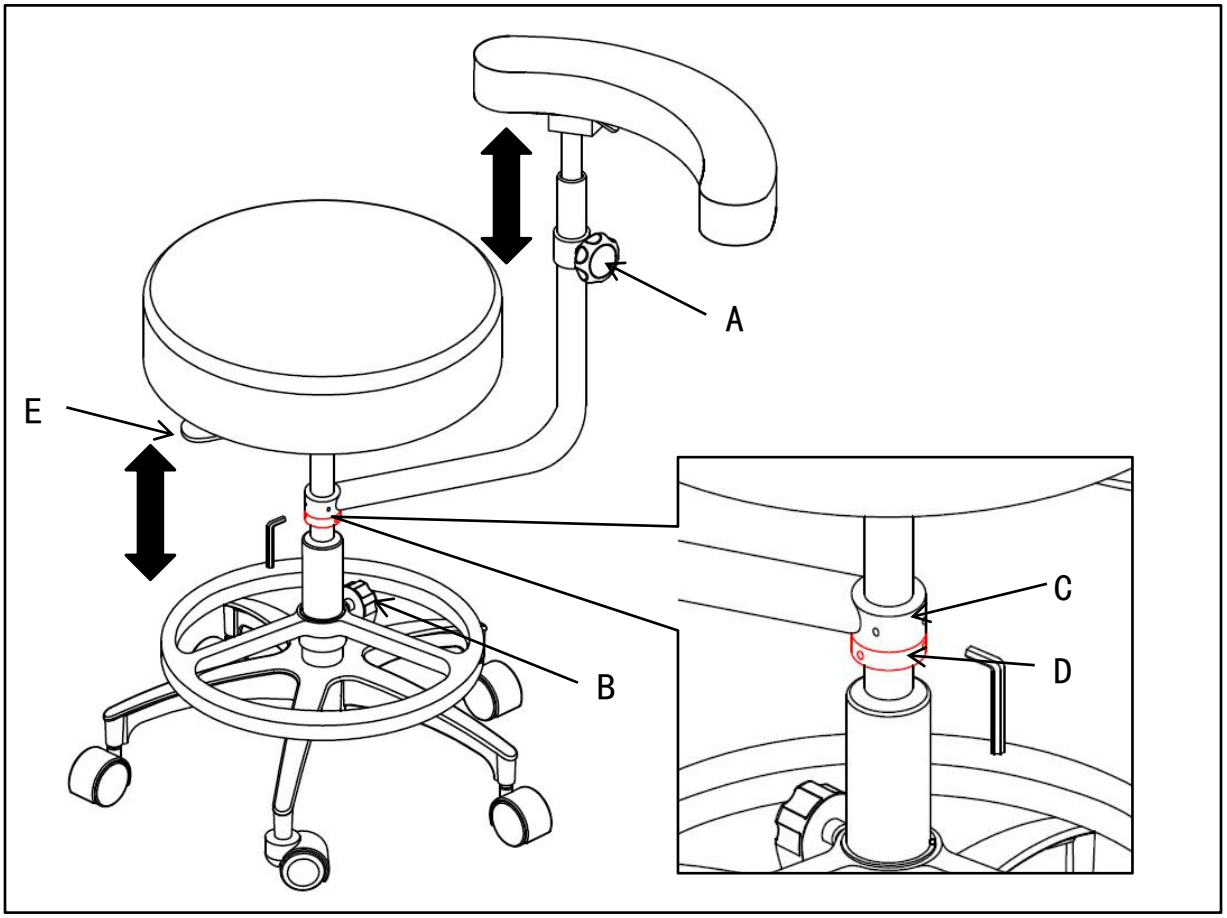


A. An adjustment handwheel of armrest.
Rotate the handwheel counterclockwise to loosen the armrest. Move the armrest to the expected position. Rotate the handwheel clockwise to secure the armrest.
B.The foot ring adjusts the handwheel.
Rotate the handwheel counterclockwise to loosen the foot ring. Move the foot ring to the expected position. Rotate the handwheel clockwise to secure the foot ring.
C. The rotation tension of the mounting arm upright can be adjusted by adjusting the screw tightness.
D. (N4 dental stool only) Height limit stop of the mounting arm upright.
Loosen the limit screw. Adjust the limit stop to the expected position. Lock the limit screw.
E. Seat height adjustment lever.
Sit on the seat. Pull up and hold the lever. Adjust the seat to the expected position. Release the lever.

3.2 Level the Armrest

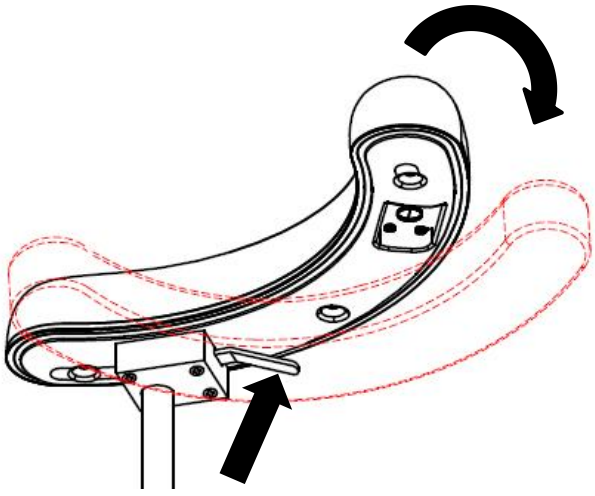

To rotate the armrest outwardly:
1. Flip and hold the armrest horizonal position lock outwardly.
2. Rotate the armrest outwardly to expected position.
3. Release the lock.

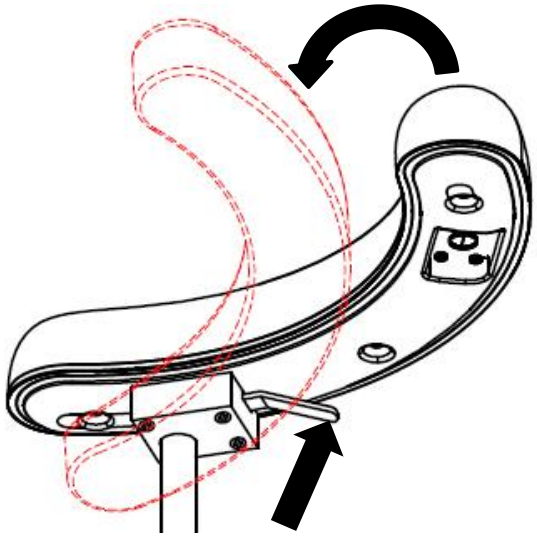

To rotate the armrest inwardly:
1. Flip and hold the armrest horizonal position lock inwardly.
2. Rotate the armrest inwardly to expected position.
3. Release the lock.

We reserve the right to make any alteration which may be due to technical improvements.





ADS Dental System Inc.

Add: 1590 S Milliken Ave., unit A, Ontario, CA 91761
Tel: (626) 620 0456
E-mail: ads@adsequip.com
http://www.adsequip.com/